

How to Use Xylazine Test Strips

1 Prepare

- Take a 5-10 mg scoop of the substance
- Add substance to about 1 teaspoon of clean water, or enough to fill the cap of a water bottle
- Mix the substance and water in a dry container (a cooker, bottle cap, or baggie the substance came in)


 ² 10 mg = enough to cover the hair of Abraham Lincoln's head on a penny

2 Test substance


- Hold the **solid red** end of the test strip
- Dip the test strip into the water up to the red line for 15 seconds then take it out
- Wait for 1-2 minutes to check your results



3 Read results

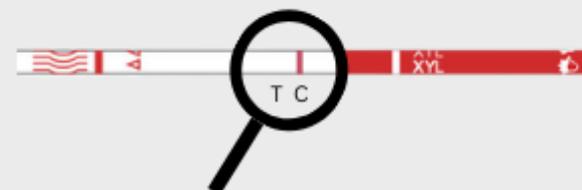

ONE red line:
Positive for xylazine
This means xylazine is present

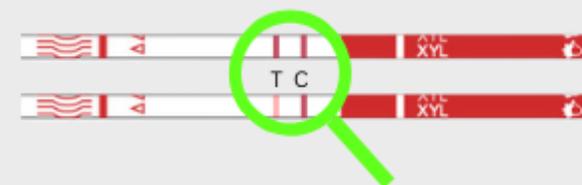

TWO red lines:
Negative for xylazine
This means xylazine is not present

The second line can be faint,
this still means negative



Scan here to visit our website for
more resources and information



MATTERS
Medication for Addiction Treatment & Electronic Referrals